

PARTS LIST Model 30911:



Enjoy the Freedom of Independant Living

- No Monitoring Fees, Service Charges or Contracts Ever
- Talk Directly through Your Pendant to a 911 Emergency Operator at the Push of a Button
- Complete Coverage In and Around Your Home
- The World's Only 2-Way Emergency Pendant Communicator
- 24 Hour 7 Day a Week Emergency Protection

LogicMark, LLC
10106 Bluegrass Pkwy
Louisville KY, 40299

One Year Limited Warranty
US Patent # 5,521,582
FCC and DOC approved
Complies with NENA guidelines
(911 Association)

1-703-934-7934
1-800-519-2419

Brochure #5001 v2.4
© 2012 LogicMark, LLC

SDG

Quick Start User's Guide



**GUARDIAN
ALERT 911**
LogicMark


GUARDIAN ALERT 911™



LogicMark

Quick Start Guide

STEP 1: Plug phone line and power cord into the back of the base.

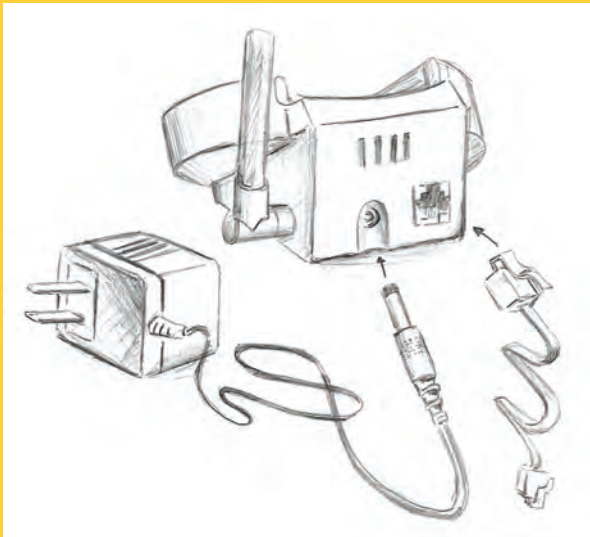


If you have any questions regarding this product, please call 1-800-519-2419 or 1-703-934-7934

STEP 2: Plug phone line and power cord into wall jacks.



If you need to plug a phone and the Guardian Alert into a single phone jack, use the splitter as shown below:

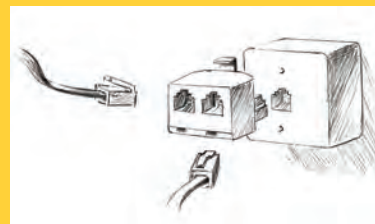

STEP 3: IMPORTANT - run a full system check on the pendant by pressing and holding the white button on the bottom for 5 seconds. You should hear a dial tone and the announcement "all systems are ok."